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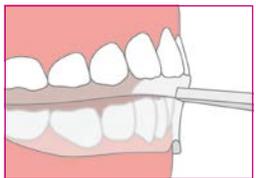
How is BRUXI+ worn? BRUXI+ – a lower jaw silicone splint is worn only at nighttime and protects tooth enamel. Oral vestibule and is automatically held in place by the lip muscles.

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Instructions for use BRUXI+ can be shortened by the dentist if necessary. Clean BRUXI+ thoroughly before using it for the first time under warm running water. Check BRUXI+ regularly for places that are bitten through! If necessary, replace BRUXI+.

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Cleaning instructions Clean it thoroughly every day with a toothbrush under warm running water. Alternately, you can also use the Dr. Hinz Orofan® dental cleaning tabs. BRUXI+ is neither dishwasher- nor microwave-safe. Clean the hygienic BRUXI+ only with a damp cloth.

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Immediate help for 3- to 12-year-old children in case of bruxism
Prevents nightly bruxism, relaxes the masticatory muscles
and protects the tooth enamel.

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Teeth grinding in children *Approximately 20 % of the children and adolescents are affected by nighttime bruxism. Noisy bruxing can already occur after eruption of the primary molars at the age of two. Causes presumably are psychological reasons such as stress, anxiety, aggressiveness or an unfavorable occlusion that needs grinding-in.*

Applied at an early stage – possibly from the age of three – BRUXI+ can effectively stop bruxing when worn at nighttime for a few weeks. Also applicable for all other age groups.

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BRUXI+ rules

Your support and a few rules will help your child:

- 1 **BRUXI+** must be voluntarily worn.
- 2 Children who have not used a pacifier must first get used to having an extraneous body in the mouth. The “chewing gum effect” of the silicone base material helps in this.
- 3 For easier habituation, **BRUXI+** should be worn in the first week one hour before going to bed and at night. If **BRUXI+** does not fall out of the mouth, nocturnal use is sufficient. If the silicone splint detaches from the teeth, wearing time should be extended (during the day by the hour while playing or watching TV).
- 4 As a rule, after wearing it regularly for about 8 -12 weeks, **BRUXI+** can be taken off on a trial basis and the grinding behavior observed. If the damaging habit continues, wearing time is extended.
- 5 After successful treatment, **BRUXI+** should nonetheless still be worn for 6 to 8 additional weeks (“danger of relapse”).
- 6 **BRUXI+** can be stored in the hygiene box.
- 7 Don't forget: Brush your teeth and **BRUXI+** on a regular basis!

If you have questions, please contact your dentist who will support the therapy!