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MUPPY® designs MUPPY® small cap is available in two different sizes:

- with red ring for deciduous dentition (size I)
- with blue ring for mixed dentition (size II)

There are also different colours and hardness grades (elastic, rigid) available.

MUPPY® Käppchen

elastisch (nur in transparent)
elastic (only transparent)

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Instructions for use Before using MUPPY® for the first time, clean it thoroughly under warm running water. Check MUPPY® regularly for places that have been bitten through! Replace MUPPY® when necessary.

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Cleaning instructions Clean thoroughly with a toothbrush under warm running water. Alternatively, you can also use the Dr. Hinz Orofan® dental cleaning sticks. MUPPY® is neither dishwasher- nor microwave-safe. Clean the hygienic MUPPY® Box only with a damp cloth.

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For straightening protruded incisors or when the
lower jaw is retracted (class II)

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MUPPY® for straightening *MUPPY® with small (biting) cap prevents tipping of the vestibular plate when the lower jaw is retracted, thereby guiding it forwards.*

Apart from the negative aesthetic impression, these malformations also have functional and health consequences: Your child cannot bite off and swallow correctly, he/she also "lisps" and breathes through the mouth instead of through the nose, which can lead to frequent colds, among other things. The goal of therapy is the weaning of habits and the mouth closure.

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The wearing calendar can be downloaded from dr-hinz-dental.de/muppy

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
MUPPY® rules *Studies have shown that more than 90 % of children give up damaging sucking with the help of MUPPY®. Your support and a few rules help your child:*

- 1 *MUPPY® is a substitute for thumb and dummy. The substitution must be voluntary.*
- 2 *MUPPY® must always be worn at night.*
- 3 *To get used to it easier children should wear MUPPY® also a few hours during the day (when playing or watching television).*
- 4 *As a rule, weaning is successful within a week. However, MUPPY® should also be worn for 6–8 weeks after successful weaning of the dummy ("recidivism").*
- 5 *MUPPY® can be stored in the hygiene box.*
- 6 *Do not forget to regularly brush the teeth and MUPPY®!*

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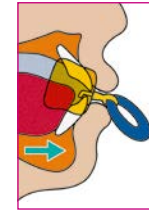How is MUPPY® worn? *MUPPY® is worn similarly to a common weaning dummy. The shield lies between the lips and the small cap behind the lower front teeth. The small cap guides the lower jaw to a neutral bite position, thereby facilitating and training the lip closure.*

If your child is having speech therapy treatment, you should take note of the following:

- Follow your speech therapist's individual wearing and exercise instructions exactly.
- Support your child during this stage, as he/she has to feel the importance of mouth closure in order to take the exercises seriously.
- Make your child aware of the "rules of play" and supervise the necessary exercises.